

## Medical Imaging Project

### Members of team:

Adriel Arias – 1621983 (soft number), Computer Science, Senior undergraduate

Jose Alvarez – 1587271 (soft number), Computer Science, Senior undergraduate

### Introduction:

The purpose of our project is to be able simulate data acquisition. By analyzing this data and applying certain algorithms to them like Fourier transforms we should be able to more accurately depict the data in a manner that is user friendly.

### Aims:

1. To develop a user friendly, yet visually pleasing, GUI.
2. Generate Phantoms for testing and validation.
3. To apply Fourier transforms and evaluate what types of images are created from Cartesian and Radial trajectories.
4. Study the effects or the reconstructed images depending on the acquisition parameters.
5. Compare the signal intensities of the input image and the final result


### Methods

To develop our GUI (1) we used PyQt5. This allowed us to add the option for the user to select what acquisition trajectory they wanted to use. Through the GUI, the user can either choose one of our generated phantoms or input their own image with a “browse” button. Furthermore, the user can see the generated k-space before and after an acquisition trajectory was applied (3), as well as how changing their parameters, with a slider, affects the outcome (4). In order to generate the phantoms, we used Skimage and created the sizes of the objects relative to that of the image with the functions `phantom1()` and `phantom2()` (2). We then used the function `shiftedDFT()` which used numpy fft functions to calculate the Fourier Transform of the image and shift it before returning. From this we can output to the GUI how the k-space looks (3). Then, we generate a mask, depending on the user trajectory selection, to multiply with the k-space and output back to the GUI (3). From this we can now generate the original image after it has been modified with the `revertDFT()` function. Finally, we are able to compare the new signal intensity with the previous signal intensity from the input image. The way


this is done is by capturing the intensities at certain pixels of the image from 0-255 depending on how intense the objects are, similar to a histogram (5)

### Results and Discussion:

(Aim 1 and Aim 2) We developed an easy to use GUI that gives the user the option to input their custom image or generate a phantom1/phantom 2 image. The user will also be given the option to choose the phantom size (512x512/1024x1024/2048x2048/4096x4096) and the acquisition trajectory. This will then show the user the selected image as a small preview.

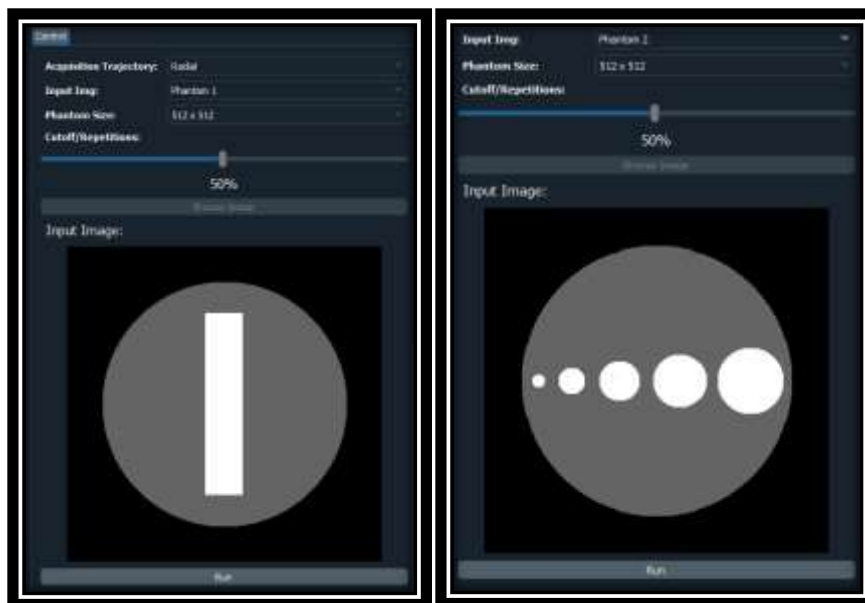


(Aim 5) Once the user presses the run button, the signal intensity of the inputted image will be presented underneath. Similarly, the signal intensity of the resulting image will be displayed at the top right of the application.

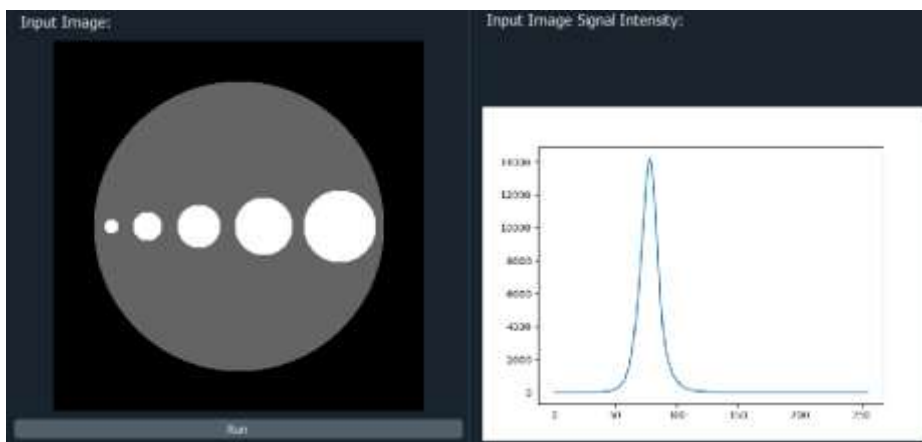


(Aim 3) In our GUI we can also display to the user the k-space generated from the input image and the collected k-space according to the acquisition trajectory they previously set. The image below shows the radial acquisition trajectory being applied with the radius being 50% of the image




This is how the final GUI interface would look to the user after the final image is calculated from the acquisition trajectory:

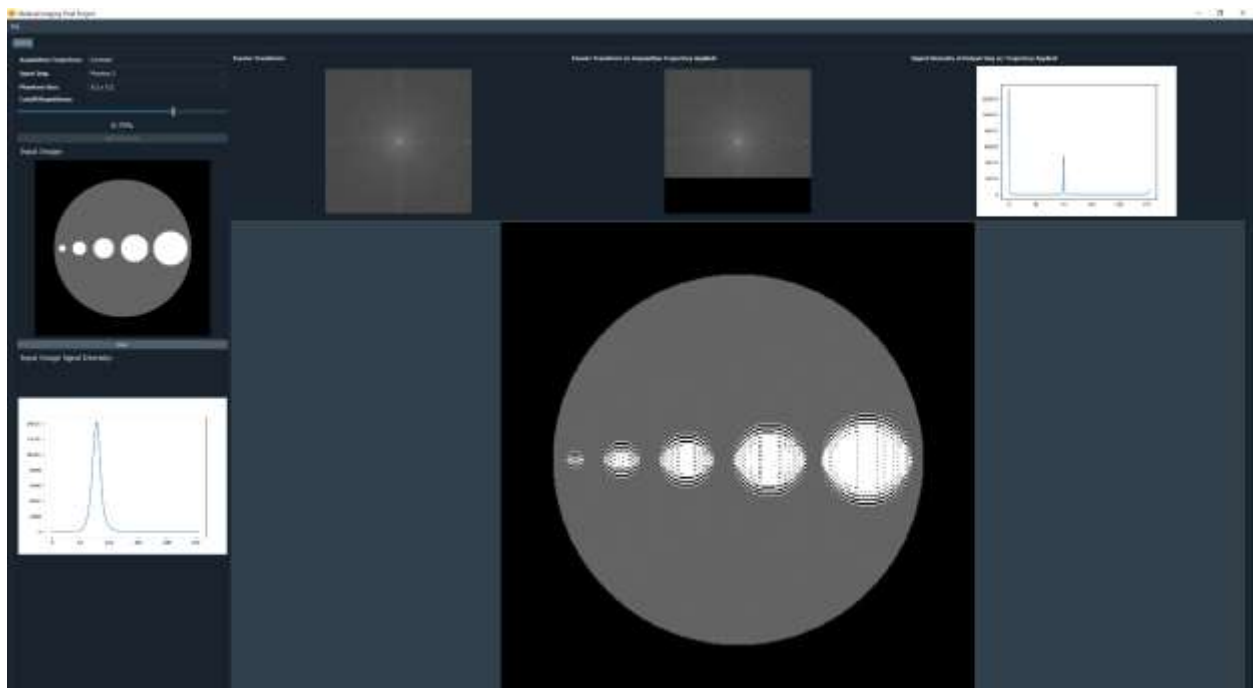


This allows the user to compare the final and original phantom. If one wants to see how changing the acquisition trajectory, one can simply change the slide under “cutoff/repetitions”. The images below demonstrate how changing this slider can impact the final result:

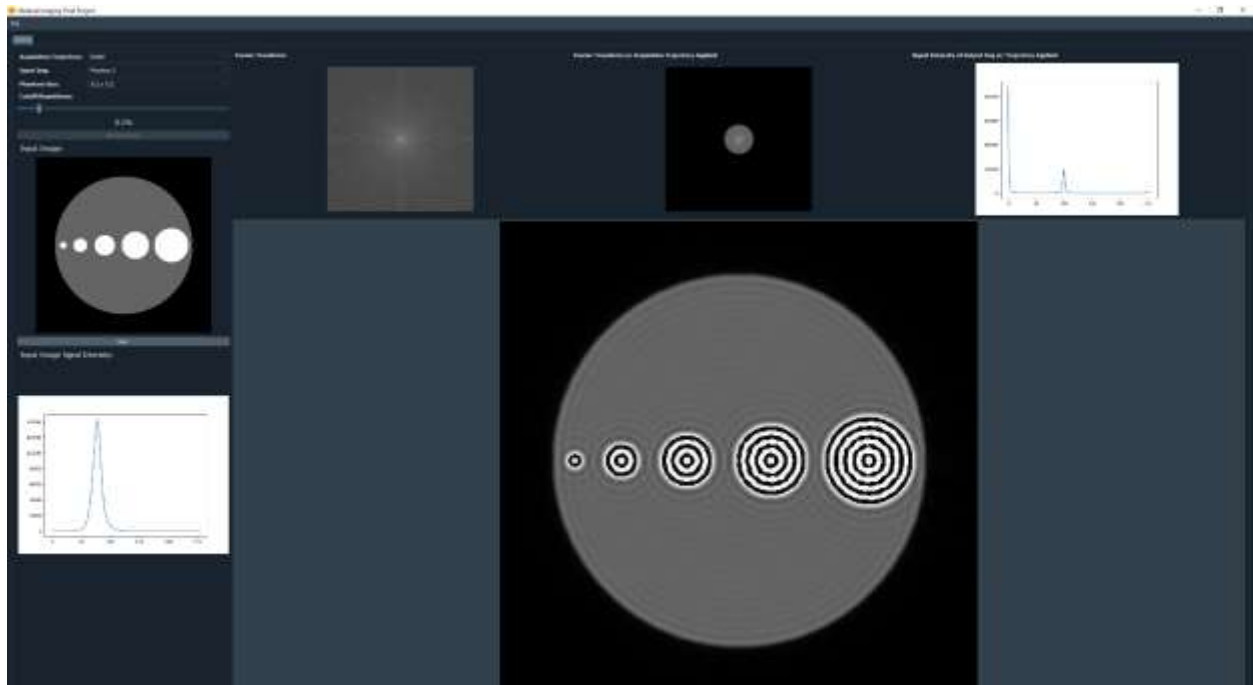

*Radial Acquisition Trajectory at 10%*

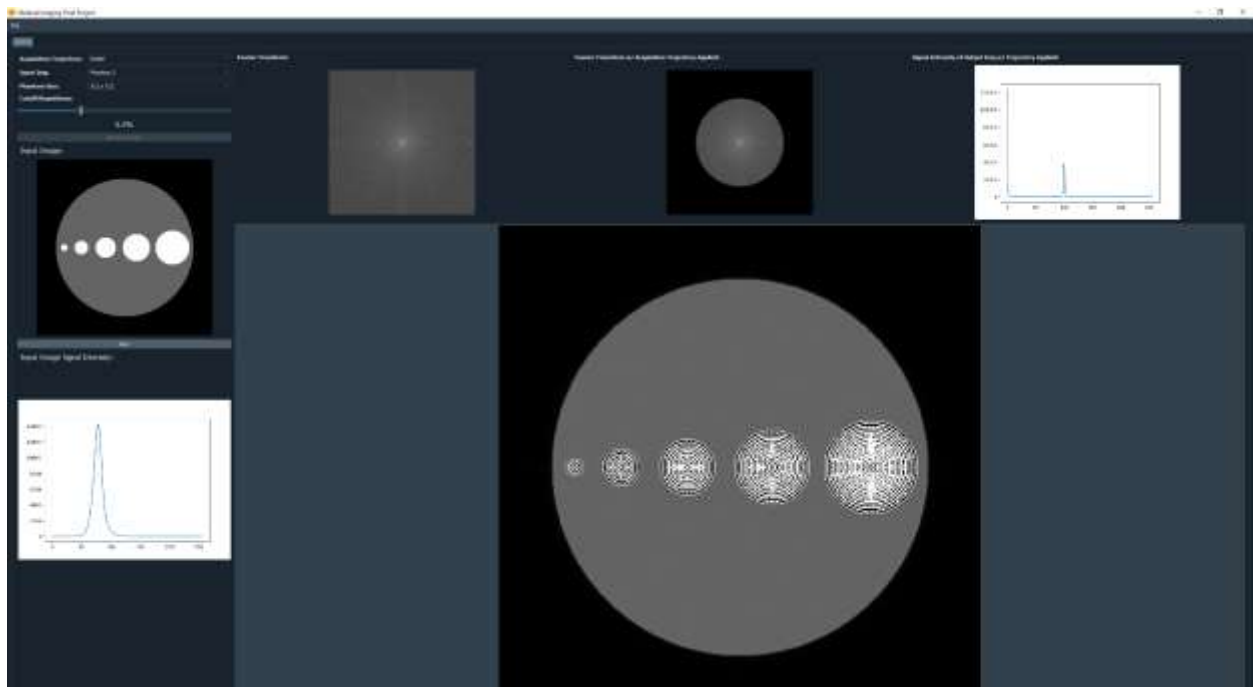


*Radial Acquisition Trajectory at 30%*

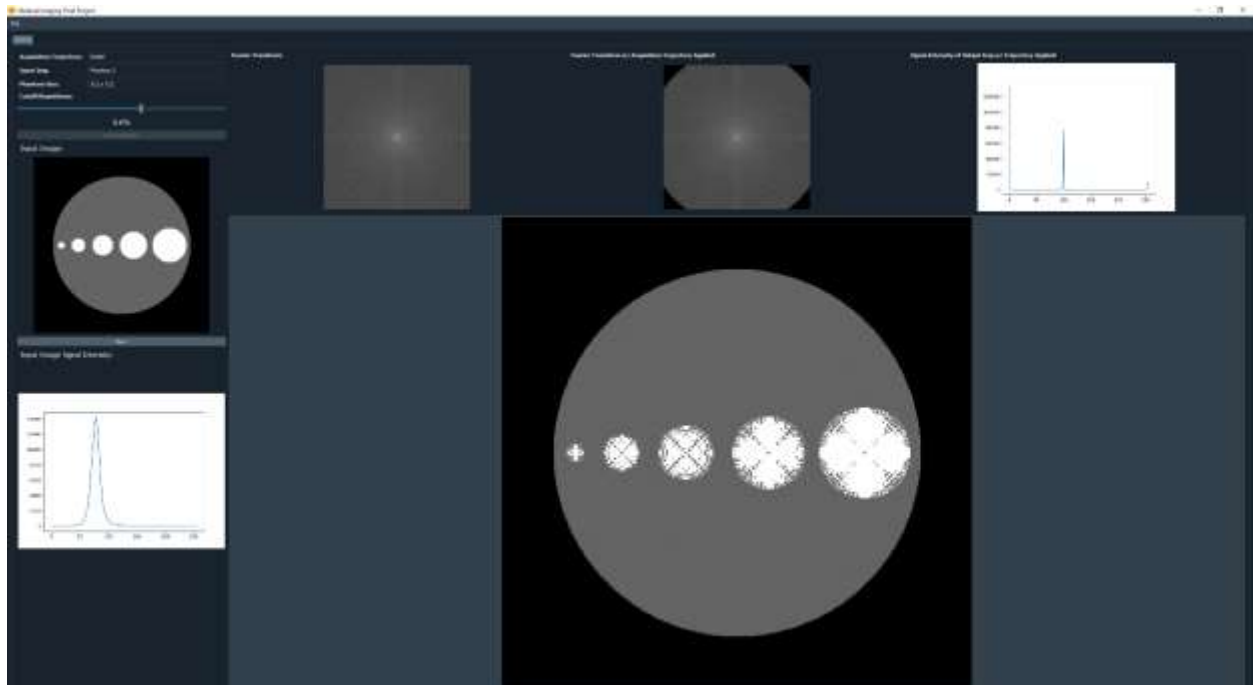

*Radial Acquisition Trajectory at 60%*

The same can be done for Cartesian!

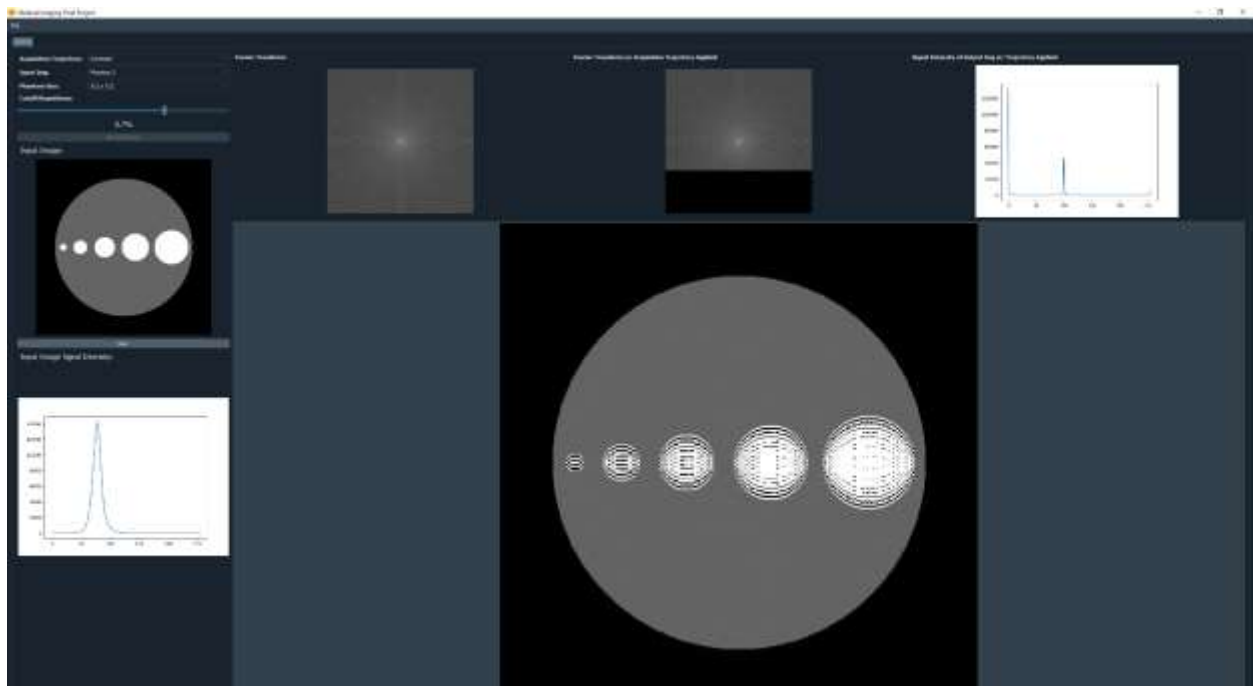

*Cartesian Acquisition Trajectory at 70%*

### Specifics:

8. Using the test/validation phantom investigate the effect of changing acquisition parameters. In all cases the "ground truth" is your original phantom! You will compare the images generated by your code relative to this phantom! Use the different tools in section 5 to investigate how your scanner works. Compare the image and the validation phantom for different acquisition parameters:
  - How the edges of the rectangular structure the first phantom are appear in the image? Any improvements?
  - In the phantom with the circular structures, what is the effect on them?
  - What is the overall image impression?
  - What parameters give you the best possible image?
  - Add a time parameter and (e.g. this will be the duration of data acquisition per horizontal line for the Cartesian acquisition, and per radial line for the radial acquisition) and investigate how "the longer the acquisition time the better the collected images"

### Answers:

- There seems to be no improvement from collecting certain k-spaces.
- The phantom with the circular structures seems to improve the longer the acquisition time is
- The overall image impression is that the parameters depend a lot on the image quality
- The parameters that give the best result are the radial acquisitions with a longer time.
- It seems that the longer the acquisition the better because it allows more data in

### **Conclusions**

We were able to successfully execute all our aims. The user is can choose what trajectory they want with the image of their choice to analyze the reconstructed images. To further analyze the reconstructed images, the user can also change the parameters and see how changing the trajectory affects the final result.

**Github link:** <https://github.com/surrealvarez/Medical-Imaging-Final>

The instructions to run the GUI can be found in the readme.

#### **CODE:**

```
import numpy as np
import cv2
from skimage.draw import circle
from skimage.draw import rectangle
from math import floor
from PIL import Image

class DFT:

    #generates the first phantom image
    def phantom1(self, x, y):
```

```

#begins the first image with the size (x,y)
img = np.zeros((y,x),dtype='uint8')

rr, cc = circle(y/2,x/2,floor(x*.39))
img[rr,cc] = 100

recHeight = floor(y*.58)
recWidth = floor(x*.12)

originY = floor((y/2)-(recHeight/2))
originX = floor((x/2)-(recWidth/2))
#coordinate of the origin
start = (originY, originX)
extent = (floor(y*.58),floor(x*.12))
rr, cc = rectangle(start, extent=extent,shape=img.shape)

img[rr,cc] = 255

# cv2.imwrite('phantom1.jpg',img)
# cv2.imread('phantom1.jpg')
return img

#generates the second phantom image
def phantom2(self, x, y):

    img = np.zeros((y,x),dtype='uint8')

    #main circle holding others
    rr, cc = circle(x*.5,y*.5,int(round(x*.390625)))

    img[rr,cc] = 100
    #generates multiple circles at a ratio

    #vertical position, horizontal position, diameter
    rr, cc = circle(int(y*.5), int(round(x*.15625)),
int(round(x*.01953125)) )
    img[rr,cc] = 255

    rr, cc = circle(int(y*.5), int(round(x*.25390625)),
int(round(x*.0390625)))
    img[rr,cc] = 255

    rr, cc = circle(int(y*.5), int(round(x*.390625)),
int(round(x*.05859375)))
    img[rr,cc] = 255

    rr, cc = circle(int(y*.5), int(round(x*.56640625)),
int(round(x*.078125)))
    img[rr,cc] = 255

    rr, cc = circle(int(y*.5), int(round(x*.771484375)),
int(round(x*.09765625)))
    img[rr,cc] = 255

    return img

```



```

def load_display(self, image):

    #Create a output window
    cv2.namedWindow("Image", cv2.WINDOW_AUTOSIZE)

    #shows image
    cv2.imshow("Image", image)

    cv2.waitKey(0)

    #destroy windows
    cv2.destroyAllWindows()

def magnitude(self, image):
    #gets the width and height of the input image
    height, width = image.shape
    #iterates and ensures values aren't negative
    for h in range(height):
        for w in range(width):
            # if negative convert to positive
            if image[h][w] < 0:
                image[h][w] *= -1
    return image

def radial(self, image, cutoff):

    """Computes a Ideal low pass mask
    takes as input:
    shape: the shape of the mask to be generated
    cutoff: the cutoff frequency of the ideal filter
    returns a ideal low pass mask"""
    height = image.shape[0] # P
    width = image.shape[1] # Q

    mask = np.copy(image)

    for u in range(height):
        for v in range(width):
            #used the pythagorean theorem to check if the coordinate
            position is within range
            duv = ((u - height/2)**2 + (v - width/2)**2)**(1/2)
            if duv <= cutoff:
                mask[u, v] = 1 #If in range leave at 1
            elif duv > cutoff:
                mask[u, v] = 0 #If out of range change to black

    maskedimg = np.multiply(image, mask)

    return maskedimg

def cartesianMask (self, image, repetitions):

    height,width = image.shape
    mask = image.copy()
    #0 is black

```

```

#iterates through image depending on the user input
for h in range (height):
    for w in range(width):
        if h < repetitions:
            mask[h,w] = 1
        elif h >= repetitions:
            mask[h,w] = 0

maskedimg = np.multiply(image, mask)

return maskedimg
def shiftedDFT(self, img):
    # run dft on image
    img = np.fft.fft2(img)
    img = np.fft.fftshift(img)
    return img

def revertDFT(self, img):
    #revert DFT to get back calculated image
    img = np.fft.ifftshift(img)
    img = np.fft.ifft2(img)
    return img

def prepareFinalOutput(self, img):
    # Steps before outputting!-----

    # this image will be returned
    magDFT = img.copy()
    magDFT = self.magnitude(magDFT)
    # magDFT = 10 * np.log(magDFT)
    magDFT = magDFT.astype("uint8") #Makes it a type that can output
    # -----
    return magDFT

def prepareOutput(self, img):
    #Output for DFT k-spaces
    # this image will be returned
    magDFT = img.copy()
    magDFT = self.magnitude(magDFT)
    magDFT = 10 * np.log(magDFT)
    magDFT = magDFT.astype("uint8") #Makes it a type that can output

    # -----
    return magDFT

```

## Code for GUI:

```

import locale
import sys
import numpy
import cv2

import qdarkstyle
from PyQt5.QtWidgets import QMainWindow
from fbs_runtime.application_context.PyQt5 import ApplicationContext

```

```

from PyQt5.QtWidgets import QApplication, QWidget, QMainWindow, QPushButton,
QAction, QMessageBox
from PyQt5.QtWidgets import QCalendarWidget, QFontDialog, QColorDialog,
QTextEdit, QFileDialog
from PyQt5.QtWidgets import QCheckBox, QProgressBar, QComboBox, QLabel,
QStyleFactory, QLineEdit, QInputDialog
from PyQt5.QtGui import QPixmap
from PyQt5 import QtGui
from window import Ui_MainWindow
from FT import DFT
from message import Ui_mssgWidget
from message import Ui_mssgWidget
from matplotlib import pyplot as plt
from PyQt5 import QtCore
# from QtCore import Qt
# from PIL import Image
import os
from scipy import stats

locale.setlocale(locale.LC_ALL, '')
# import operator
path = os.path.dirname(os.path.abspath(__file__))

class Main(QMainWindow, Ui_MainWindow):

    def __init__(self):
        super(Main, self).__init__()
        self.setupUi(self)

        self.actionQuit.triggered.connect(QApplication.quit)
        self.browseButton.clicked.connect(self.openFile) # This is used as a
signal to browse for an img
        self.runButton.clicked.connect(self.onClickRun) # This runs the
function onClickRun when run button is clicked
        self.functionSelected = str(self.comboBoxTrajectory.currentText()) #
Current trajectory selected
        self.imgSelect = str(self.comboBoxInputImg.currentText()) # Current
img Choice selected
        self.phantomSizeSelected = str(self.comboBoxSize.currentText()) #
Select size of Phantom
        # print(self.functionSelected)
        # self.comboBox.activated[str].connect(self.onComboClicked)
        #
self.comboBoxHistShapingParams.activated[str].connect(self.onComboClicked2)
        self.comboBoxInputImg.currentTextChanged.connect(self.onComboClicked)
        #
self.histScaleSlider.valueChanged.connect(self.histScaleValue.setNum)
        self.cutoffSlider.valueChanged.connect(self.cutoffLabel)
        self.cutoffValue = float(self.cutoffSlider.value())/100

        self.MRI = DFT()

    def compute_histogram(self, image):

        hist = [0] * 256

```

```

r, c = image.shape

for i in range(r):
    for j in range(c):
        hist[image[i, j]] = hist[image[i, j]] + 1

return hist

def cutoffLabel(self,value):
    self.cutoffValue = float(value)/100
    self.cutoff_Ratio.setText(str(self.cutoffValue) + '%')

# # -----
# ##### THIS METHOD IS IN CHARGE OF KEEPING TRACK OF WHICH FUNCTION IS
SELECTED IN THE DROP DOWN MENY #####
# def onComboClicked2(self, text):
#     # Here, we change the value of global
#     self.functionSelectedParam =
str(self.comboBoxHistShapingParams.currentText())
#     print(self.functionSelectedParam)
#
#####
#####
#
# # -----
# ##### THIS METHOD IS IN CHARGE OF KEEPING TRACK OF WHICH FUNCTION IS
SELECTED IN THE DROP DOWN MENU #####
def onComboClicked(self):
    self.inPic.clear() # We first clear the img
    self.name = '' # RESETS self.name
    inputDir = path + '/Inputs/'

    if not os.path.exists(inputDir):
        os.makedirs(inputDir)

    # Here, we change the value of global
    self.inputImgSelection = str(self.comboBoxInputImg.currentText())
    # IF a phantom is selected, we disable the browse button
    if self.inputImgSelection == 'Phantom 1':
        self.name = '' # RESETS self.name
        inputDirP1 = inputDir + 'phantom1.jpg'
        self.browseButton.setEnabled(False)
        img = self.MRI.phantom1(360,360)
        cv2.imwrite(os.path.join(inputDir,'phantom1.jpg'), img)

        self.pixmap = QPixmap(appctxt.get_resource(inputDirP1))
        self.inPic.setPixmap(self.pixmap.scaled(self.inPic.width(),
self.inPic.height(), QtCore.Qt.KeepAspectRatio))
        self.inPic.setAlignment(QtCore.Qt.AlignCenter)

    elif self.inputImgSelection == 'Phantom 2':
        self.name = '' # RESETS self.name
        inputDirP2 = inputDir + 'phantom2.jpg'
        self.browseButton.setEnabled(False)
        img = self.MRI.phantom2(360,360)

```

```

        cv2.imwrite(os.path.join(inputDir,'phantom2.jpg'), img)

        self.pixmap = QPixmap(appctxt.get_resource(inputDirP2))
        self.inPic.setPixmap(self.pixmap.scaled(self.inPic.width(),
self.inPic.height(), QtCore.Qt.KeepAspectRatio))
        self.inPic.setAlignment(QtCore.Qt.AlignCenter)
        # ELSE, the browse button is available
    else:
        self.browseButton.setEnabled(True)

    #
#####
#####

    # -----
-----

    ##### THIS FUNCTION BROWSES AND SAVES THE PIXMAP OF INPUT IMG AND
    FILEPATH OF INPUT IMAGE AS self.name #####
    def openFile(self):

        try:
            self.name, _ = QFileDialog.getOpenFileName(self, 'Open File',
options=QFileDialog.DontUseNativeDialog)
            print(self.name, ' Input Img')
            self.pixmap = QPixmap(appctxt.get_resource(self.name))
            self.inPic.setPixmap(self.pixmap.scaled(self.inPic.width(),
self.inPic.height(), QtCore.Qt.KeepAspectRatio))
            self.inPic.setAlignment(QtCore.Qt.AlignCenter)

        except:
            pass

#####
#####

    def onClickRun(self):

        cutoffFactor = float(self.cutoffValue)
        outputDir = 'Output/'
        outputDir2 = path + '/Output/'

        if not os.path.exists(outputDir2):
            os.makedirs(outputDir2)

        # Here we do the operation for the Cartesian function

        try:
            self.inputHist.clear()
            self.resultHist.clear()
            self.inputImgSI.clear()
            self.outImgSI.clear()
        except:
            print("ERROR!!!")

        # If the chosen img is Phantom1, we then define the size by the
        selected Size
        if str(self.comboBoxInputImg.currentText()) == 'Phantom 1':

```

```

        if(str(self.comboBoxSize.currentText()) == '512 x 512'):
            # print(512)
            size = 512
            img = self.MRI.phantom1(512,512)
        elif(str(self.comboBoxSize.currentText()) == '1028 x 1028'):
            # print(1028)
            size = 1028
            img = self.MRI.phantom1(1028, 1028)
        elif(str(self.comboBoxSize.currentText()) == '2056 x 2056'):
            # print(2056)
            size = 2056
            img = self.MRI.phantom1(2056, 2056)
        elif(str(self.comboBoxSize.currentText()) == '4096 x 4096'):
            # print(4096)
            size = 4096
            img = self.MRI.phantom1(4096,4096)
        # If the chosen img is Phantom1, we then define the size by the
        selected Size
        elif str(self.comboBoxInputImg.currentText()) == 'Phantom 2':
            if(str(self.comboBoxSize.currentText()) == '512 x 512'):
                # print(512)
                size = 512
                img = self.MRI.phantom2(512,512)
            elif(str(self.comboBoxSize.currentText()) == '1028 x 1028'):
                # print(1028)
                size = 1028
                img = self.MRI.phantom2(1028, 1028)
            elif(str(self.comboBoxSize.currentText()) == '2056 x 2056'):
                # print(2056)
                size = 2056
                img = self.MRI.phantom2(2056, 2056)
            elif(str(self.comboBoxSize.currentText()) == '4096 x 4096
(Computationally Heavy)'):
                # print(4096)
                size = 4096
                img = self.MRI.phantom2(4096,4096)
        # If the chosen imgs is Custom, we browse for an img
        elif str(self.comboBoxInputImg.currentText()) == "Custom Img (Use
Browse)":
            try:
                img = cv2.imread(self.name, 0)
                height, width = img.shape
                size = height
            except:
                self.mssgBox = QWidget()
                message = Ui_mssgWidget()
                message.setupUi(self.mssgBox)
                message.mssgLabel.setText("Browse for an\image first!")
                self.mssgBox.setWindowFlags(QtCore.Qt.Window |
QtCore.Qt.CustomizeWindowHint |
QtCore.Qt.WindowTitleHint |
QtCore.Qt.WindowCloseButtonHint | QtCore.Qt.WindowStaysOnTopHint)
                self.mssgBox.show()
            return

        # height, width = img.shape

```

```

        img = self.MRI.shiftedDFT(img)
        kSpaceImgOut = self.MRI.prepareOutput(img) # Fourier transform Output
with Cutoff
        cv2.imwrite(os.path.join(outputDir2, 'kSpaceImg.jpg'), kSpaceImgOut)

        kSpaceImgOut = cv2.imread(outputDir2 + '/kSpaceImg.jpg', 0)
        inHist = self.compute_histogram(kSpaceImgOut)
        # TO COMPUTE SI AND DISPLAY
        plt.clf()
        inHist_fig = plt.plot(inHist)
        plt.savefig(outputDir2 + '/Input_SI_Graph.png')
        self.pixmapInSI = QPixmap(outputDir2 + '/Input_SI_Graph.png')

self.inputImgSI.setPixmap(self.pixmapInSI.scaled(self.inputImgSI.width(),
self.inputImgSI.height(), QtCore.Qt.KeepAspectRatio))
        self.inputImgSI.setAlignment(QtCore.Qt.AlignCenter)

        self.pixmapInHist = QPixmap(outputDir2 + '/kSpaceImg.jpg')

self.inputHist.setPixmap(self.pixmapInHist.scaled(self.inputHist.width(),
self.inputHist.height(), QtCore.Qt.KeepAspectRatio))
        self.inputHist.setAlignment(QtCore.Qt.AlignCenter)

        cutoff = cutoffFactor * size

        if str(self.comboBoxTrajectory.currentText()) == "Cartesian":
            outCartesianPath = path + '/Output/Cartesian_Output/'
            if not os.path.exists(outCartesianPath):
                os.makedirs(outCartesianPath)

                kSpaceImgCut = self.MRI.cartesianMask(img, cutoff)
                kSpaceImgOutCut = self.MRI.prepareOutput(kSpaceImgCut) # Fourier
transform Output with Cutoff

cv2.imwrite(os.path.join(outCartesianPath, 'kSpaceImgCut.jpg'), kSpaceImgOutCut
)
            # kSpaceImgOutCut = cv2.imread(outCartesianPath +
'/kSpaceImgCut.jpg', 0)
            self.pixmapResultHist = QPixmap(outCartesianPath +
'/kSpaceImgCut.jpg')

self.resultHist.setPixmap(self.pixmapResultHist.scaled(self.resultHist.width(
), self.resultHist.height(), QtCore.Qt.KeepAspectRatio))
        self.resultHist.setAlignment(QtCore.Qt.AlignCenter)

        outImgCut = self.MRI.revertDFT(kSpaceImgCut)

        outImgCut = self.MRI.prepareFinalOutput(outImgCut) #

cv2.imwrite(os.path.join(outCartesianPath, 'kSpaceImgCut.jpg'), outImgCut)

        outImgCut = cv2.imread(outCartesianPath + 'kSpaceImgCut.jpg', 0)
        plt.clf()
        outHist = self.compute_histogram(outImgCut)
        inHist_fig = plt.plot(outHist)
        plt.savefig(outCartesianPath + '/Output_SI_Graph.png')

```

```

        self.pixmapOutSI = QPixmap(outCartesianPath +
'/Output_SI_Graph.png')

self.outImgSI.setPixmap(self.pixmapOutSI.scaled(self.outImgSI.width(),
self.outImgSI.height(), QtCore.Qt.KeepAspectRatio))
self.outImgSI.setAlignment(QtCore.Qt.AlignCenter)

self.pixmapOut = QPixmap(outCartesianPath + 'kSpaceImgCut.jpg')
# self.outPic.setScaledContents(True)
self.outPic.setPixmap(self.pixmapOut.scaled(self.outPic.width(),
self.outPic.height(), QtCore.Qt.KeepAspectRatio))
self.outPic.setAlignment(QtCore.Qt.AlignCenter)

elif str(self.comboBoxTrajectory.currentText()) == "Radial":
    outRadialPath = path + '/Output/Radial_Output/'
    if not os.path.exists(outRadialPath):
        os.makedirs(outRadialPath)

    kSpaceImgCut = self.MRI.radial(img, cutoff)
    kSpaceImgOutCut = self.MRI.prepareOutput(kSpaceImgCut) # Fourier
transform Output with Cutoff

cv2.imwrite(os.path.join(outRadialPath, 'kSpaceImgCut.jpg'), kSpaceImgOutCut)
# kSpaceImgOutCut = cv2.imread(outRadialPath +
'/kSpaceImgCut.jpg', 0)
self.pixmapResultHist = QPixmap(outRadialPath +
'/kSpaceImgCut.jpg')

self.resultHist.setPixmap(self.pixmapResultHist.scaled(self.resultHist.width(
), self.resultHist.height(), QtCore.Qt.KeepAspectRatio))
self.resultHist.setAlignment(QtCore.Qt.AlignCenter)

outImgCut = self.MRI.revertDFT(kSpaceImgCut)
outImgCut = self.MRI.prepareFinalOutput(outImgCut) #

cv2.imwrite(os.path.join(outRadialPath, 'kSpaceImgCut.jpg'), outImgCut)
outImgCut = cv2.imread(outRadialPath + 'kSpaceImgCut.jpg', 0)
outHist = self.compute_histogram(outImgCut)
plt.clf()
inHist_fig = plt.plot(outHist)
plt.savefig(outRadialPath + '/Output_SI_Graph.png')
self.pixmapOutSI = QPixmap(outRadialPath +
'/Output_SI_Graph.png')

self.outImgSI.setPixmap(self.pixmapOutSI.scaled(self.outImgSI.width(),
self.outImgSI.height(), QtCore.Qt.KeepAspectRatio))
self.outImgSI.setAlignment(QtCore.Qt.AlignCenter)

self.pixmapOut = QPixmap(outRadialPath + 'kSpaceImgCut.jpg')
# self.outPic.setScaledContents(True)
self.outPic.setPixmap(self.pixmapOut.scaled(self.outPic.width(),
self.outPic.height(), QtCore.Qt.KeepAspectRatio))
self.outPic.setAlignment(QtCore.Qt.AlignCenter)

```



```
if __name__ == '__main__':  
    appctxt = ApplicationContext()           # 1. Instantiate ApplicationContext  
    window = Main()  
    window.showMaximized()  
    appctxt.app.setStyleSheet(qdarkstyle.load_stylesheet_pyqt5())  
    exit_code = appctxt.app.exec_()          # 2. Invoke appctxt.app.exec_()  
    sys.exit(exit_code)
```